


Popular (/category/337/snl-10)

SNL-10

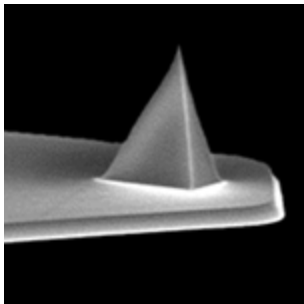

A

B

C

D

Geometry		Tip Radius (nm)	
▼ B Triangular		Nom: 2	
		Max: 12	
Frequency (KHz)		Length (µm)	
Nom: 23		Nom: 205	
Min: 16		Min: 200	
Max: 28		Max: 210	
Spring Const (N/m)		Width (µm)	
Nom: 0.12		Nom: 40	
Min: 0.06		Min: 35	
Max: 0.24		Max: 45	

Related Products



Order a free SNL-10 sample

Order a Free Sample



(/) Please note that we are only able to offer samples for shipping addresses within the United States of America.

\$324.50 (USD)

* Sold in packs of 10

1

Add to Cart

Buy more & save!



SHIPPING INFO (T-SHIPPING.ASPX)



PROBE SELECTION GUIDE (T-APPLICATION-SELECTOR-GUIDE.ASPX)



EMAIL US (MAILTO:AFMPROBEORDERS@BRUKER.COM)

Questions? Free, Online Consulting **+1 (800) 715-8440 (tel:+1 (800) 715-8440)**

Overview

Bruker's Sharp Nitride Lever (SNL) probes are the industry standard for high-performance imaging with TappingMode in fluid. These probes marry the sharpness of a Silicon tip with the low spring constant and high sensitivity of a Silicon Nitride cantilever, for an unprecedented level of resolution and force control. The unique design of each SNL probe includes four cantilevers with varying degrees of force constants and resonant frequencies, providing enough versatility to cover a large majority of sample types.

The SNL cantilever layout is identical to that of the other NP products, with "A" and "B" cantilevers on one side of the probe, and "C" and "D" cantilevers on the other side of the probe.

See the cantilever orientation diagram [here](#).

(<https://www.brukerafmprobes.com/Images/SiN%20Sheet.jpg>) All cantilevers on the SNL


BRUKER

probes have ± 2 degrees of cantilever bend, making them ideal for use on all Bruker and non-Bruker AFMs.



(/)

Tip Specification

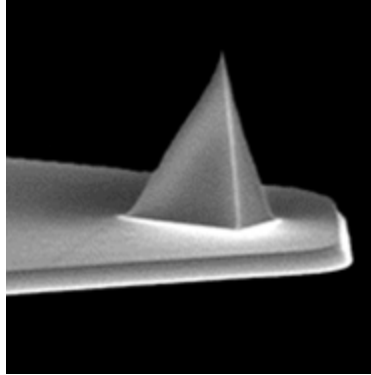

 [View Larger](#)


Tip Height (h): 2.5 - 8.0 μm

Front Angle (FA): $15 \pm 2.5^\circ$

Back Angle (BA): $25 \pm 2.5^\circ$

Side Angle (SA): $22.5 \pm 2.5^\circ$

Cantilever Specification

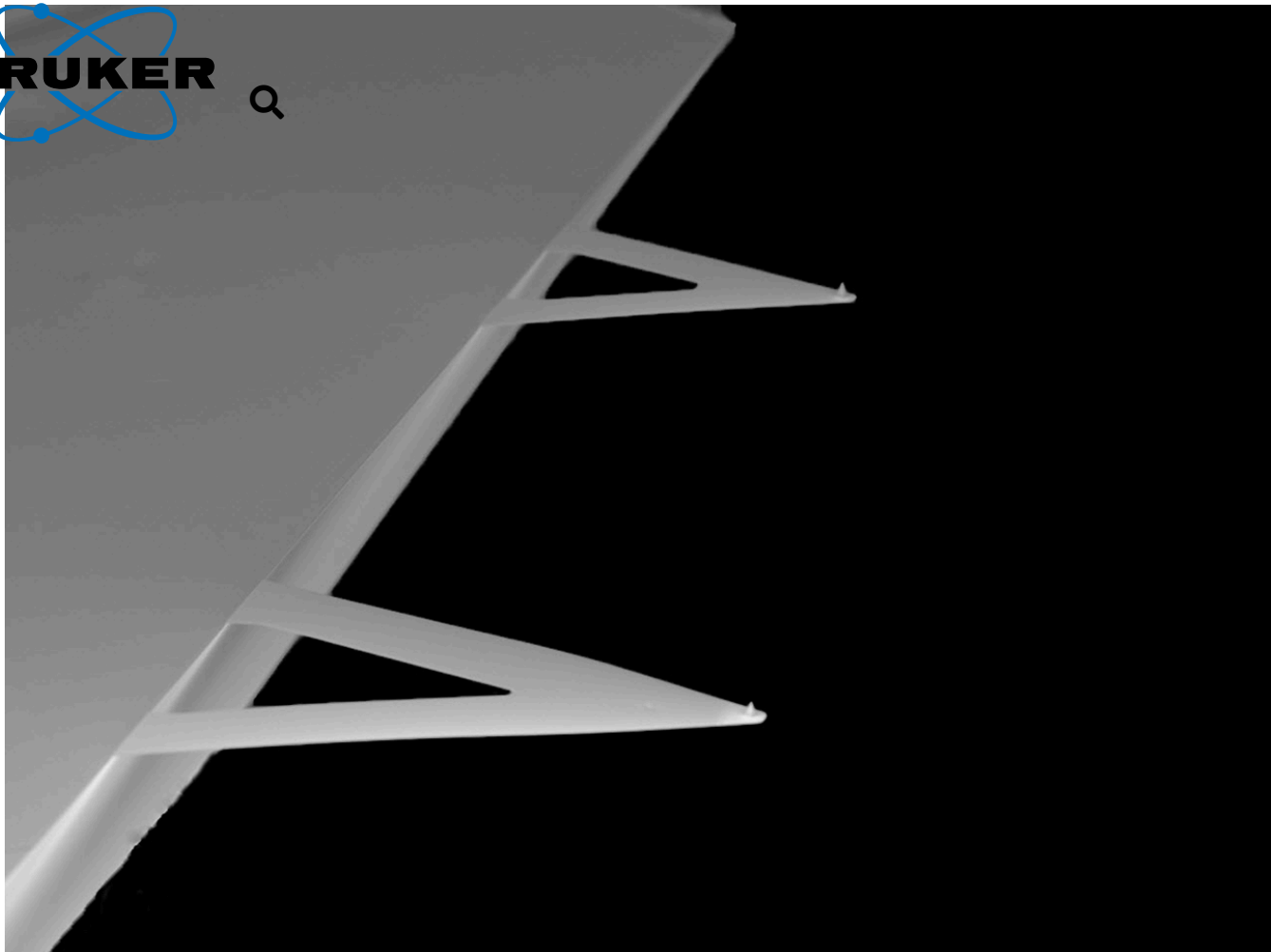

 [View Larger](#)  [Cantilever Schematic](#)

Material: Silicon Nitride

Geometry: Triangular

Cantilevers Number: 4

Cantilever Thickness (Nom): 0.6 μ m

Cantilever Thickness (RNG): 0.55 - 0.65 μ m

Back Side Coating: Reflective Gold

Top Layer Back: 45 $\pm$ 5 nm of Ti/Au